


NOTES: UNLESS OTHERWISE SPECIFIED


1. APPLICABLE STANDARDS AND SPECIFICATIONS; ASME Y14.5-2009, DIMENSIONING AND TOLERANCING ASME Y14.5-2007, ABBREVIATIONS

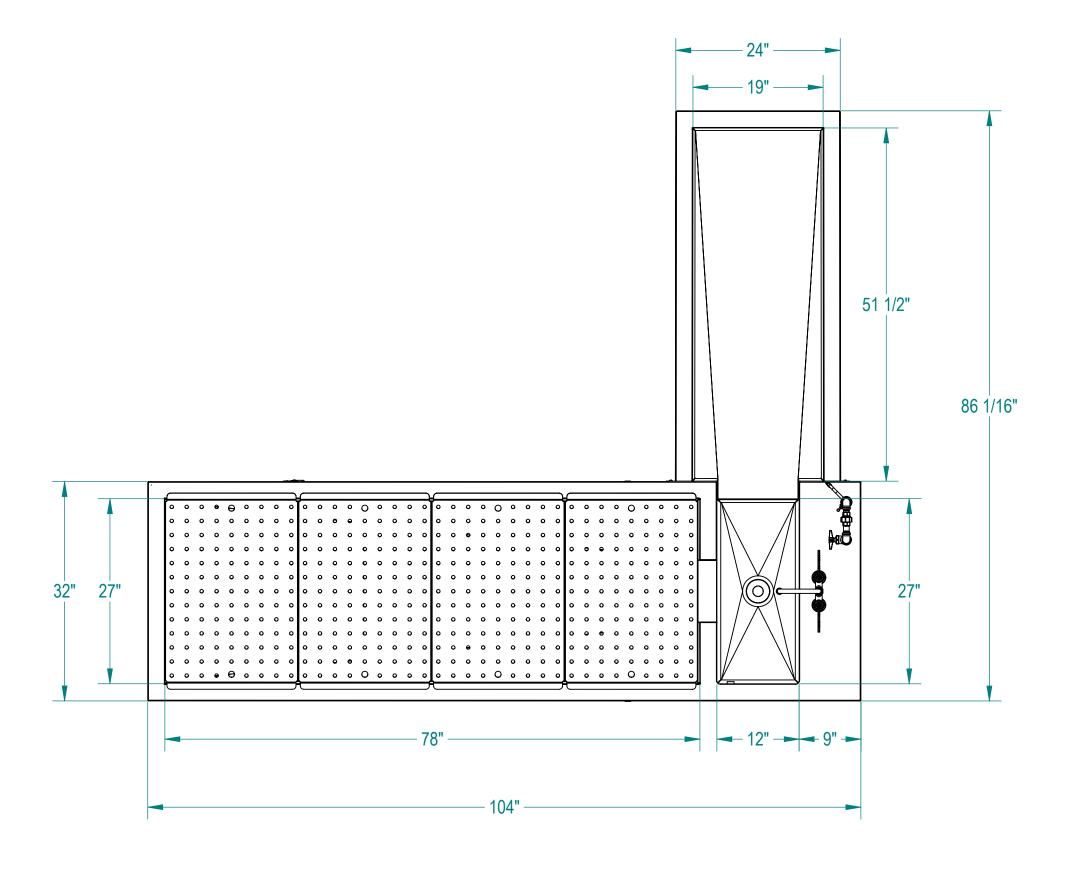

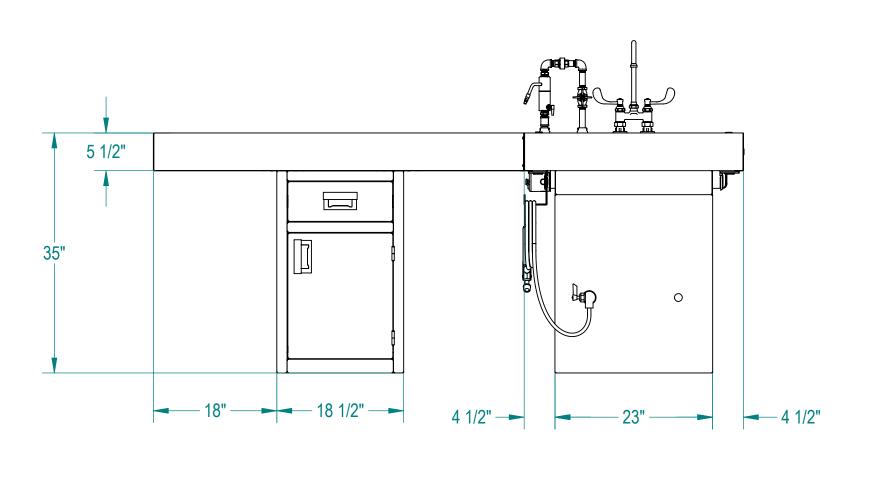

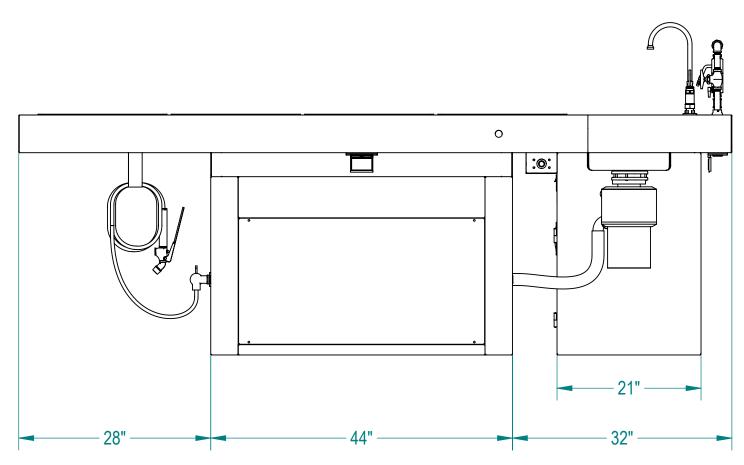


	CUSTOMER	APPROVAL		PROPRIETARY THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN			MOPEC GROUP	ELEVATING PATHOLOGY MORTECH ADVISION OF MORE CRICIL	411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com	
	PROJECT			AUTHORIZATION FROM MOPEC GROUP.			MODEL 1036-19L (CE700)			
		OLOMATURE.	DATE	DO NOT SCALE DRAWING	APPROVALS	DATE	DEDECTAL ALITO	DOV TABLE I CHADED I	I CHYDED LEET	
				MATERIAL	DRAWN BY	11/05/04	PEDESTAL AUTOPSY TABLE, L-SHAPED, LEFT			
		SIGNATURE	DATE	304 STAINLESS STEEL	ELISEO S. CHECKED BY	11/05/24	SIZE DWG NO.	1036-19L	REV.	
				FINISH						
				#4	APPROVED BY		SCALE 1:14 NEXT ASSY.		SHEET 1 OF 4	

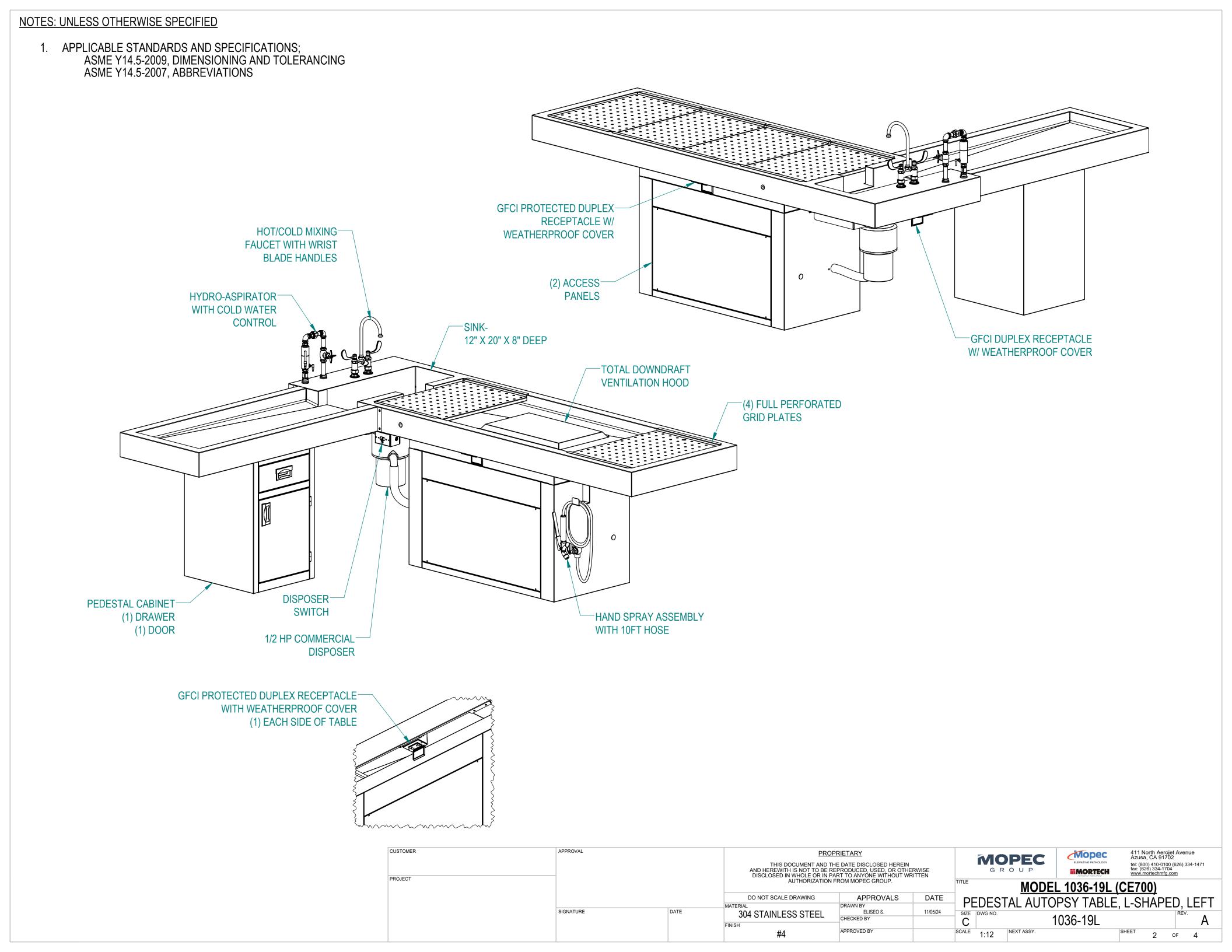

UTILITY REQUIREMENTS:


- 1. A 4" X 5" AREA FOR INLET PLUMBING TO INCLUDE THE FOLLOWING:
 - (1) COLD WATER SUPPLY, 3/4" OD COPPER STUBBED OUT 6"
 - (1) HOT WATER SUPPLY, 1/2" OD COPPER STUBBED OUT 6"
 - SUPPLY SHUT-OFF VALVES AT EACH WATER SUPPLY FLUSH LINES CLEAR BEFORE CONNECTING TO STATION
- 2. (1) STANDARD 1 1/2" WASTE LINE STUBBED OUT 6" SUPPLY P-TRAP FOR DRAIN
- 3. A 4" X 4" AREA FOR ELECTRICAL INPUT TO INCLUDE THE FOLLOWING:
 - (1) 115V/60Hz/1PH/15A CIRCUIT FOR STATION
 - (1)115V/60Hz/1PH/15A CIRCUIT FOR DISPOSER
- 4. A 6" DIA DUCT EXHAUST CONNECTION:
 - RECOMMENDED EXHAUST: 450-500 CFM AT 0.50 STATIC PRESSURE MEASURED AT THE 6" EXHAUST, PROVIDE WATER CAP OVER DUCT INLET.
 - EXISTING CONDITIONS MAY WARRANT AN INCREASE OR DECREASE, CONTACT YOUR ENGINEERING DEPT.
- 5. PROVIDE 60-65 PSI WATER PRESSURE AND 4 GALLONS PER MINUTE OF FLOW FOR PROPER OPERATION OF HYDRO ASPIRATOR
- 6. ELECTRICIAN TO PROVIDE LIQUID TIGHT JUNCTION BOX AND CONDUIT PER LOCAL CODES.

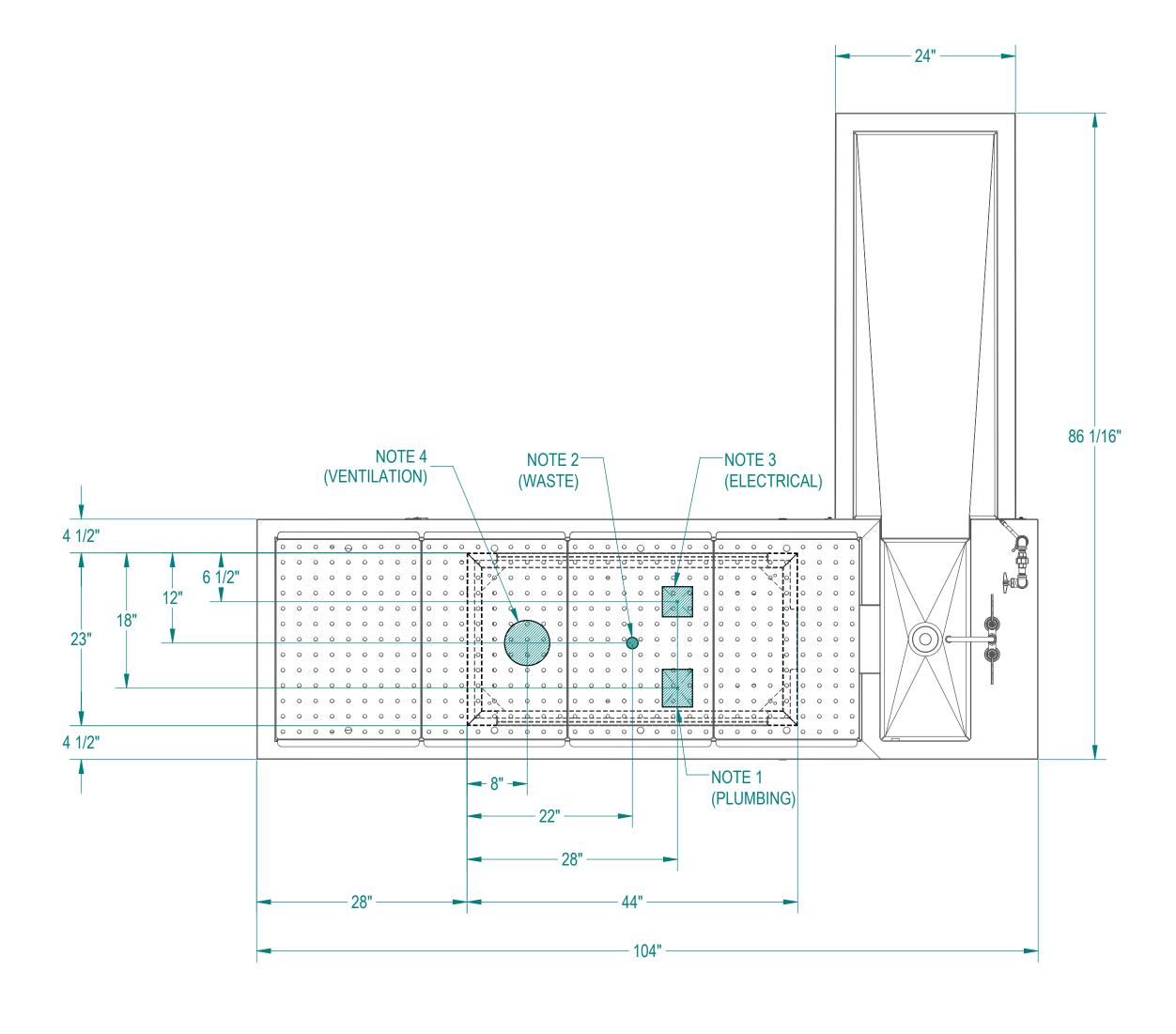

CUSTOMER	APPROVAL		THIS DOCUMENT AND TH AND HEREWITH IS NOT TO BE RI DISCLOSED IN WHOLE OR IN PA	PROPRIETARY THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MOPEC GROUP.			ELEVATING PATHOLOGY ***MORTECH A DEVICENCY IN MORE CALCUIT	411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com
PROJECT			AUTHORIZATION				MODEL 1036-19L (CE700)	
			DO NOT SCALE DRAWING	APPROVALS	DATE	PEDESTAL AUTO	PSY TARI F	, L-SHAPED, LEFT
			MATERIAL	DRAWN BY		ILDESTALAST		, L-OIIAI LD, LLI I
	SIGNATURE	DATE	304 STAINLESS STEEL	ELISEO S. CHECKED BY	11/05/24	SIZE DWG NO.	1036-19L	REV.
			FINISH				1000-13L	
			#1	APPROVED BY		SCALE 1.12 NEXT ASSY.		SHEET 3 OF 1

NOTES: UNLESS OTHERWISE SPECIFIED

1. APPLICABLE STANDARDS AND SPECIFICATIONS; ASME Y14.5-2009, DIMENSIONING AND TOLERANCING ASME Y14.5-2007, ABBREVIATIONS

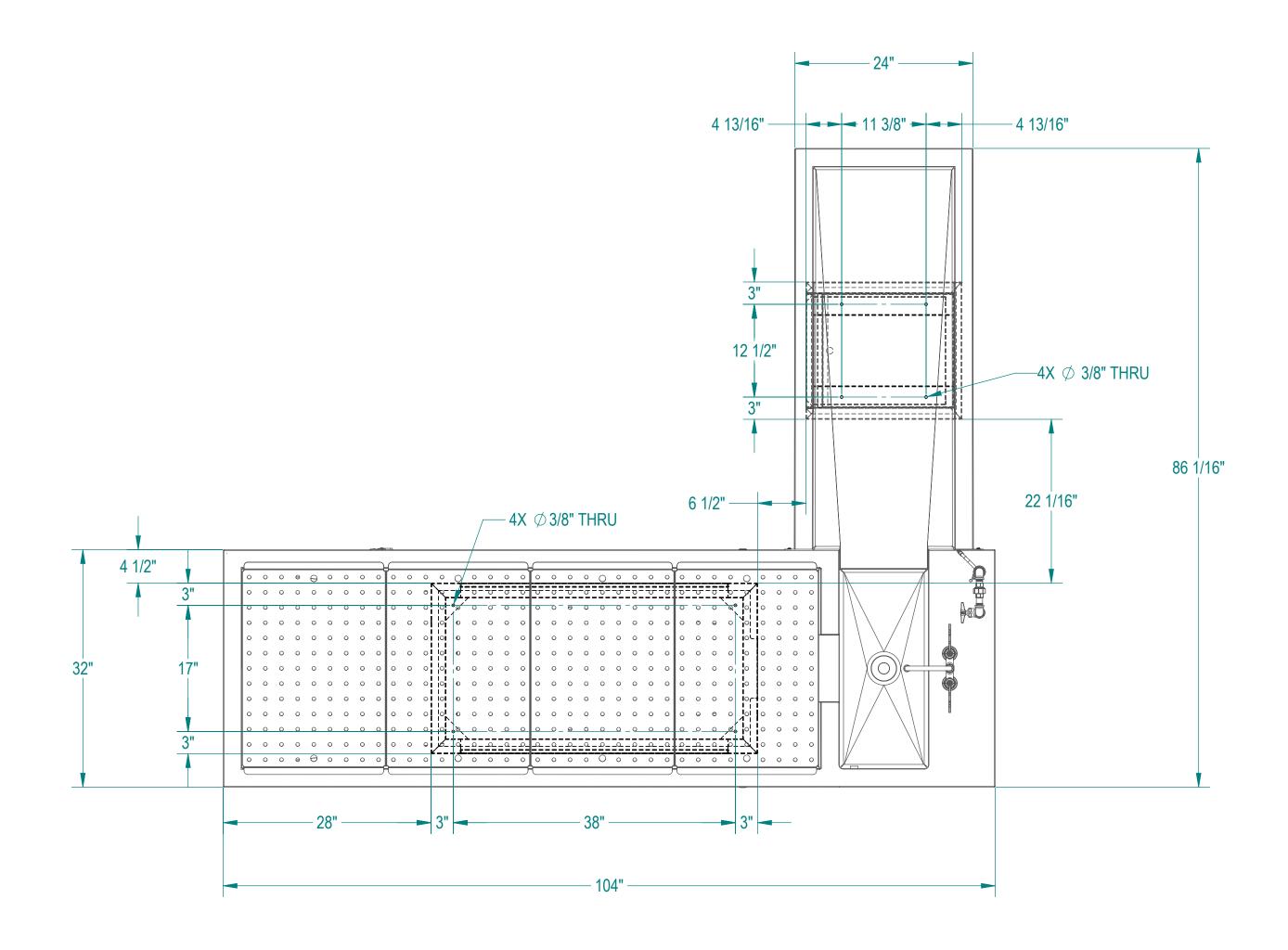

FLOOR MOUNTING DETAIL

